

Applications


Amplifier seeding
Coherent THz generation
Multiphoton microscopy
Time-resolved spectroscopy
Asynchronous optical sampling

FEMTOSOURCE™ twin synergy™



FEMTOSOURCE™ twin synergy™ is the combination of two FS synergy™ oscillators on one temperature stabilized base plate, sharing one pump source. This setup results in extremely stable repetition rates of both lasers relative to each other (*Optics Express*, Vol. 13, No. 22, pp. 9029, October 31, 2003). Each FS SYNERGY™ can be configured individually. The ultrafast Ti:Sapphire oscillators employ our patented Dispersive Mirror (DM) technology, which is specifically designed for broadband intracavity group delay dispersion compensation.

The system generates shortest optical pulses with unprecedented quality, stability, and reproducibility. DM systems constitute the ideal solution for a wide range of scientific, industrial, and medical applications.

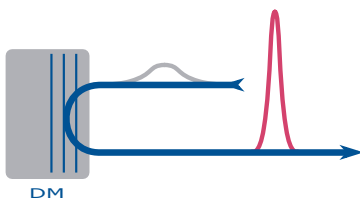

FEMTOSOURCE™

twin synergy™


Extraordinary Features

- Minimal drift between both cavities
- Extremely broad spectrum
- Extremely low noise
- Excellent long term output stability
- Ultrashort pulses | Compact footprint


Options

- Length stabilisation with FEMTOLOCK™
- Custom repetition rate
- CE-Phase stabilisation

TWIN SYNERGY™ - FRONT VIEW, Dimensions in [mm] ([in])

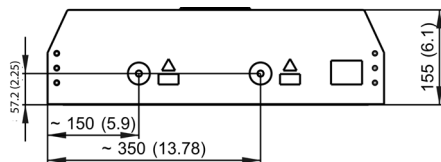

TWIN SYNERGY™ - TOP VIEW, Dimensions in [mm] ([in])

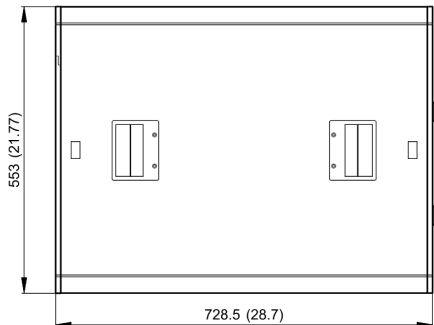

Ultrashort Pulses

Intracavity high order dispersion can be eliminated to an extent that high quality, bandwidth limited sub-10 fs pulses can be generated from a Ti:Sapphire oscillator (*Optics Letters, Vol. 20, No. 6, pp. 602, March 15, 1995*). An Extra Cavity Dispersion Control (ECDC) unit provides precise control of dispersion outside of the oscillator to transmit pulses with these extraordinary parameters to your target.

Compactness

DM systems contain no intracavity elements other than the gain medium resulting in an extremely compact and simple setup.

Stability & Reproducibility

The intracavity dispersion is not sensitive to cavity alignment, in strong contrast to prism controlled systems. Hence the day-to-day reproducibility as well as the stability of the laser output parameters is excellent in the sub-10 fs regime. Owing to the compact DM resonator design the system features the lowest energy noise and timing jitter demonstrated in the sub-10 fs regime to date.

Flexibility

Custom designed systems are available with repetition rates between 65 and 125 MHz. The center wavelength can be chosen and factory set to your specifications. Each Oscillator can be configured as FS synergy™ PRO or as FS synergy™ 20. The output power for each oscillator can be chosen at factory depending on pump power distribution.

FEMTOSOURCE™	TWIN SYNERGY™ PRO	TWIN SYNERGY™ 20
Pulse duration	< 10 fs	< 20 fs
Bandwidth (FWHM) @ 800 nm	> 100 nm	> 40 nm
Mode-locked output power (av.)	> 400 mW @ 5 W (532 nm) ²⁾	> 400 mW @ 5 W (532 nm) ²⁾
Output energy @ 75 MHz	> 5.3 nJ	> 5.3 nJ
Peak power @ 75 MHz	> 530 kW	> 265 kW
Beam diameter (1/e ²)	< 2 mm	
Spatial mode	TEM ₀₀ (M ² < 1.3)	
Polarization	> 100:1 (horiz.)	
Noise (measured 10 Hz - 100 KHz)	< 0.05 % rms	
Power stability ¹⁾	± 1%	
<small>1) Measured over 2 hours after 30 min. warmup at constant environmental conditions 2) 8 W Pump source required</small>		
<small>All specifications are subject to change without notice</small>		



VISIBLE AND INVISIBLE LASER RADIATION
 AVOID EYE OR SKIN EXPOSURE TO
 DIRECT OR SCATTERED RADIATION

LASER TYPE: Ti:Sapphire
 WAVELENGTH: 650 - 1050 nm
 PULSE DURATION: > 8 fs
 MAX. PULSE ENERGY: 20 nJ @ 75 MHz
CLASS 4 LASER PRODUCT

Synergy 1, 2008



FEMTOLASERS Produktions GmbH
 Fernkorngasse 10 | 1100 Wien | Austria
 P: +43 1 503 7002 0
 F: +43 1 503 7002 99
 info@femtolasers.com

FEMTOLASERS, Inc.
 124 Mt. Auburn St. | Suite 200 North
 Cambridge | MA 02138 | USA
 P: +1 978 456 9920
 F: +1 978 456 9922
 info@femtolasers.com



FEMTOLASERS' laser products are certified to comply with the Federal Regulations (21 CFR Subchapter J) as administered by Center of Devices and Radiological Health on all systems ordered for shipment after October 1, 2003.